


GLASS NUTSCHE FILTER

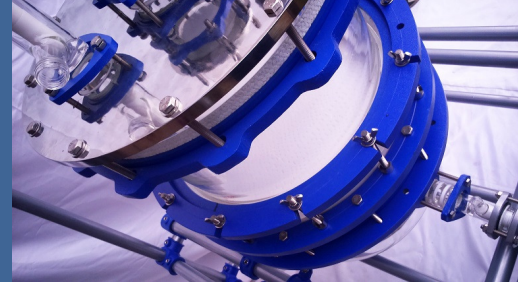

Good Partner - Good Success

❖ SALIENT FEATURES :

- ✓ Lift system for easy change of filter material.
- ✓ Tested for full vacuum and overpressure.
- ✓ Hand stirrer or optional explosion proof. electric stirrer.
- ✓ Stainless steel frame with wheels
- ✓ High Filtration Efficiency

Ablaze Nutsche Filter are specially designed for effective filtration and can be conveniently combined with our standard reaction units to enable solid-liquid separation.

All contact parts are inert and hence allow the same equipment to be used with a wide range of chemicals and solvents. Appropriate filter is used to separate the solid particles of the slurry, forming a filtered cake. The filtrate drains to the bottom and can be collected from the bottom outlet valve. The filter cake can be easily removed from the top with the help of quick release coupling.

Capacity: Up to 100Litre

Pressure: Up to +1 Bar G

Temperature: Up to 200°C

Filter Cloth: Various Mesh Size and MOC available as per requirement

Material of Construction: Borosilicate Glass, PTFE Customized Configurations available as per requirement.



Ablaze Export Incorporation
E-52/4, Sardar Estate , Ajwa Road
Vadodara 390019- Gujarat INDIA
www.ablazeexport.com